




# FACIAL PROTECTION

# THE RIGHT FACIAL PROTECTION: PUTTING YOUR SAFETY FIRST

The most important element in selecting a face mask is the confidence that comes from knowing that both you and your patients are adequately protected.

Not all face masks are created equal; some masks have a higher level of protection than others. Here at Halyard, we have displayed the American Society for Testing and Materials (ASTM) International Standards for face masks on our packaging, so that you can face the unknown, knowing you're protected.

Medical Face Mask Material Requirements by Performance LEVEL ASTM International Standard F2100-19 <sup>1</sup>			
CHARACTERISTIC	LEVEL 1 BARRIER	LEVEL 2 BARRIER	LEVEL 3 BARRIER
Bacterial filtration efficiency, %	≥ 95	≥ 98	≥ 98
Differential pressure, mm H <sub>2</sub> O/cm <sup>2</sup>	< 5.0	< 6.0	< 6.0
Sub-micron particulate filtration efficiency at 0.1 micron, %	≥ 95	≥ 98	≥ 98
Resistance to penetration by synthetic blood, minimum pressure in mm Hg for pass result	80	120	160
Flame spread	Class 1	Class 1	Class 1

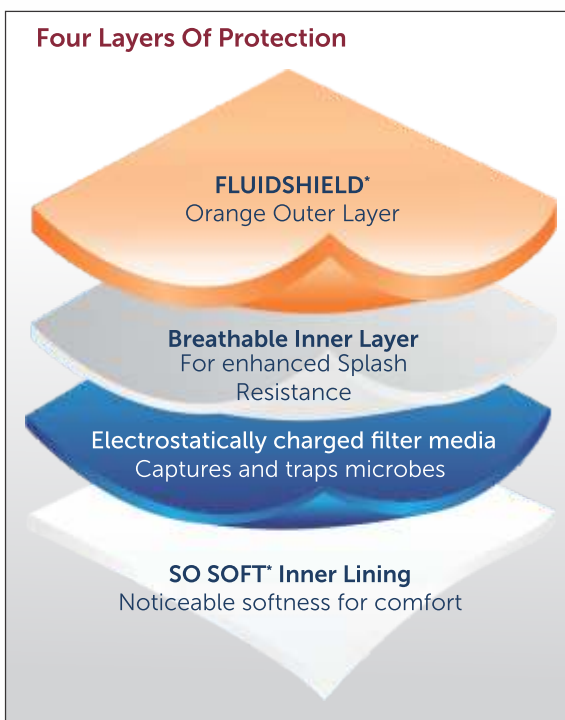

## SUPERIOR PROTECTION NEVER FELT SO GOOD

All Halyard next generation FLUIDSHIELD\* Face Masks incorporate four layers of barrier protection, including an electrostatically charged layer of filter media, for enhanced protection and an inner SO SOFT\* lining to deliver the critical balance of protection and comfort.

**You can be confident there is a Halyard\* Face Mask that meets your clinical needs and personal preferences**

1. ASTM F2100-19: <https://www.astm.org/DATABASE.CART/HISTORICAL/F2100-19.htm>



## FLUIDSHIELD\* LEVEL 3 FACE MASKS

- Fluid resistance at 160mm Hg
- Meets ASTM F2100-19 Level 3 standards for fluid resistance, filtration and breathability
- Electrostatically charged filter media for protection
- Features smooth SO SOFT\* inner lining for superior comfort



48207



48247

Item Code	Description	Quantity
48207	FLUIDSHIELD* Level 3 Fog-Free Surgical Mask with SO SOFT* Lining, Foam Band, Orange, Ties	50 masks/Box, 6 boxes/case
48247	FLUIDSHIELD* Level 3 Fog-Free Surgical Mask with SO SOFT* Lining, Wrap-around Visor, Foam Band, Ties, Orange	25 masks/Box, 4 boxes/case



## N95 PARTICULATE FILTER RESPIRATORS

- NIOSH approved N95 Particulate filter Respirator
- For procedures/situations that require Respiratory and/or Facial Splash Protection
- Fluid resistant at 160 mm Hg (Level 3), 120mm Hg (Level 2)
- Polyurethane headband
- Tested against 7 chemotherapy drugs
- Duckbill style



46727

Item Code	Description	Quantity
46767	FLUIDSHIELD* Level 3 N95 Particulate Filter Respirator and Surgical Mask with Safety Seal and SO SOFT* Lining, Regular Size, Orange	35 masks/Box, 6 boxes/case
46727	FLUIDSHIELD* Level 3 N95 Particulate filter Respirator and Surgical Mask with SO SOFT* Lining, Regular Size, Orange	35 masks/box, 6 boxes/case
62126	FLUIDSHIELD* Level 2 N95 Particulate Filter Respirator and Surgical Mask with SO SOFT* Lining, Regular Size, White	50 masks/Box 6 boxes/case

## STANDARD FACE MASKS

### Surgical Masks

- Protective three-layer construction for when exposure to blood and/or bodily fluids is not a risk
- Appropriate for general patient care and isolation



49214

Item Code	Description	Quantity
49214	Fog-Free Surgical Mask, Pleat-Style with Ties, Blue	Masks/Box = 50 Boxes/Case = 6

# FACIAL PROTECTION

## EYEWEAR

### Safeview\* Glasses

- For use in situations which may involve splashing or spraying
- Designed to be lightweight for added comfort in use
- Shipped with colorful frames to appeal to a variety of preferences
- Single-use lenses may be easily replaced into reusable frames
- Designed to be easily worn over most prescription glasses
- Constructed from natural rubber latex-free materials



<b>Code</b>	<b>Description</b>	<b>Packaging</b>
SV50A	SAFEVIEW* Assembled Glasses	5 Pkg. of 10 Each Case = 50 Assembled Glasses

## FACE SHIELD



### Medichoice Full Face Shield

- Designed for procedures and environments where risk of exposure to blood and/or bodily fluids through splashing is possible
- Equipped with a foam band and headband
- Can be worn with or without glasses or goggles
- Resistant to fogging

<b>Code</b>	<b>Description</b>	<b>Packaging</b>
FS1100F	Shield Face Full Length Visor with Anti-Fog and Anti-Glare Coating Elastic Headband Foam Header Plastic Clear	24 face shields/Box 10 boxes/case



**O&M Halyard India Pvt Ltd**  
World Trade Center, Tower-1 Unit 308  
Kharadi, Pune 411 014  
Maharashtra, India

For feedback/suggestions/complaints, contact Customer care officer  
Tel no: +91 20-46845500/5555  
Email: [customercare\\_IND@hyh.com](mailto:customercare_IND@hyh.com)

COPY-03816

O&M Halyard makes no representation that its products described in this catalogue will provide complete protection from exposure to contaminants, fluids or micro-organisms or will eliminate the risk of contracting infectious diseases.

\*Registered Trademark or Trademark of O&M Halyard Inc. or its affiliates.  
©2021. All rights reserved

[www.halyardhealth.in](http://www.halyardhealth.in)  
 O&M Halyard APAC